

Commissioning and Installation of Coordinate Detector for Super BigBite Spectrometer



Jacob Bird¹, Peter Monaghan^{1, 2}, and Ralph Marinaro^{1, 2}

¹Christopher Newport University


²Thomas Jefferson National Accelerator Facility



Abstract


Hall A at Jefferson Lab is operating a modular electromagnetic spectrometer, called Super BigBite (SBS), to investigate the charge distribution within nucleons. SBS uses a coordinate detector (CDet) to track the scattered electron position. The detector measures vertical displacement to determine the out-of-plane angle of the electron's trajectory.

Detector Design

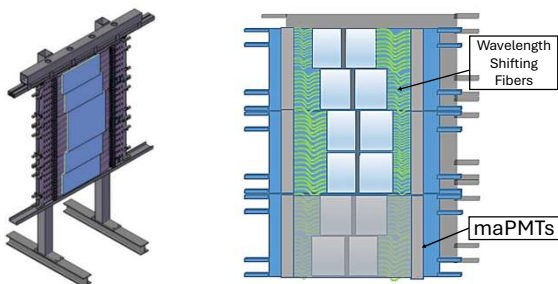

An integral aspect of the detector design is the use of Multi-Anode Photomultiplier Tubes (maPMTs). These convert photons into charge signal, which is later used for analog (ADC) and logical (TDC) signals.

The figure to the left shows the structure of one of the maPMTs. Green wavelength fibers deliver the light signal to each pixel of the maPMT, and the charge output to each corresponding channel is located on the back of the device (gray). These PMTs are powered with high voltage.

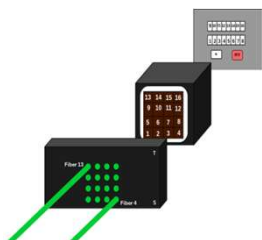


Crosstalk Methodology

A primary issue when collecting data is crosstalk between pixels in the maPMTs. A crosstalk event occurs when a signal in one channel is also registered in another channel due to capacitive coupling within the PMT. Crosstalk events are identified by their characteristically short time widths. This is observed with a multi-input oscilloscope (below).




Having recognized the short TDC widths, a cut is performed to analytically reduce the number of crosstalk events in each PMT channel (left). Fitting a plot of TDC leading vs. trailing edge helps determine the final value of the TDC width cut. This cut is applied to TDC width spectra (right) for confirmation.

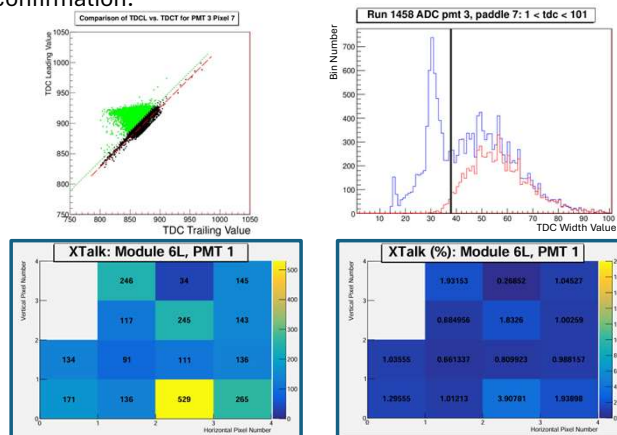

Heat maps (above) are generated to calculate how much of the collected comics data is crosstalk per pixel.

Results

Analysis was run for the first PMT on the left side of Module 6 (PMT1, M6L). On each pixel of PMT1, there are an average of 161.67 crosstalk events, with extreme values of 529 hits and 34 hits. This PMT has a total hit average per pixel of 13,446.29, meaning 1.33% of events in PMT1 are crosstalk events. This percentage closely matches the manufacturer's crosstalk value of 1%. Crosstalk analysis concluded with the generation of a plot to use when the experiment is running in Hall A. This plot gives useful information for all PMTs on Layer 1 of M6L including total hits, crosstalk hits, and double pixel hits.

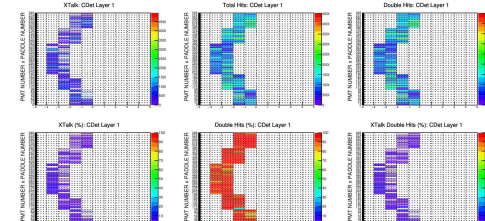

Conclusion and Future Work

The successful installation of CDet will support running experiments with SBS and further investigation of nucleon structure. Future work on this project includes final installation of CDet into Hall A, expanding all other analysis code to be detector wide, installing connectors to additional ribbon cables, and complete calibration analysis after installation.

Acknowledgements

Sincere thanks to my mentor, Dr. Peter Monaghan, and alternate mentor, Dr. Ralph Marinaro, for their support and guidance throughout this project.